Multisample Moisture Sorption Analyzer – SPS

Dynamic water vapor sorption (DVS)

Dynamische Wasserdampfsorption

Moisture

Virtually all materials interact more or less strongly with humidity. Knowledge of moisture-induced changes of the material properties is an important key parameter for decisions regarding processing, packaging, storage and shelf-life of a product.

Dynamic water vapor sorption

The gravimetric DVS method measures the change in mass over time of a sample kept in an environment of controlled and constant temperature and relative humidity. Change in mass takes place either by water uptake (sorption) from the air or water release (desorption) from the sample.

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Sorption isotherm and sorption kinetics

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The relationship between the equilibrium water content of a sample and the relative humidity of the ambient air at a particular temperature is known as the water vapor sorption isotherm.



The sorption isotherm enables to draw conclusions on fundamental material properties and provides information on material handling with respect to moisture.

From the sorption kinetics curve of a material, moisture induced structural changes are determined that are often accompanied by sudden, temporary changes in moisture content.

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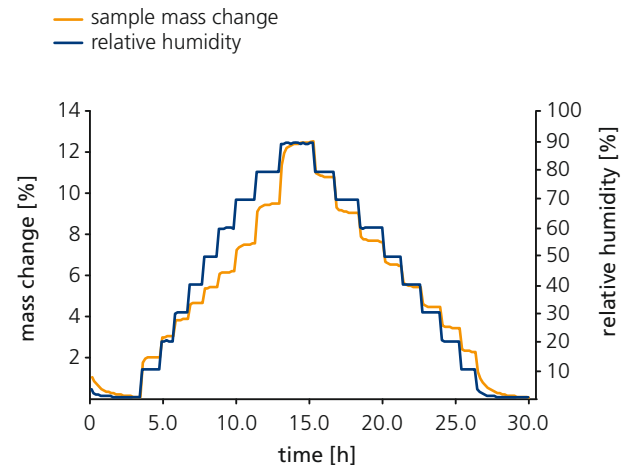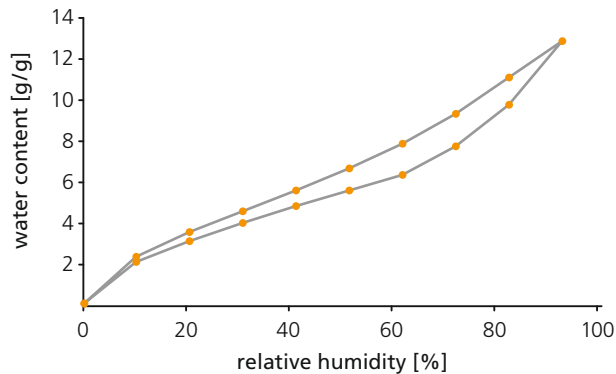Application areas

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The instruments of the SPS series differ in precision and load range of their micro- and sub-micro-balances as well as in number and size of the sample positions.

The SPS model series covers all areas of application, from pharmaceutical pre-formulation to building materials.

The choice of instrument for an application depends, among other parameters, on the size required for a representative sample, the amount of sample available and the expected total water sorption of the materials.

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SPS11-10 μ	SPSx-1 μ High Load	SPSx-1 μ Advance	SPS23-100n
			Pharm. Pre-formulation
			Pharm. Applications, (Hydrate Formation, Amorphous Content)
			Material Characterization, Polymers, Electronics
			Food Research, Food Technology, Shelf-Life
			Fine Chemicals, Detergents, Personal Care
			Natural Fibres, Textiles, Wood, Paper
			Building/Construction Materials, Fertilizer
			Permeability, Packaging

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SPS model	SPSx-1 μ Advance	SPSx-1 μ High Load
Number of samples & sample size	Exchangeable samples trays for 11 samples \varnothing 50 mm and 23 samples \varnothing 33 mm	
Temperature range	5 °C to 60 °C (\pm 0.1 K)	
Humidity range*	0 % RH to 98 % RH (\pm 0.5 % RH at 10...30 °C)	
Load range min to max	<10 mg to 22 g	<30 mg to 220 g
Balance resolution	1 μ g	1 μ g (up to 22 g) 10 μ g (up to 220 g)
Repeatability**	\pm 2 μ g	\pm 5 μ g (up to 22 g) \pm 20 μ g (up to 220 g)
Hardware options	Video camera Raman spectroscopy	Video camera Raman spectroscopy Water vapor transmission rate Large objects sample tray
Software options	21 CFR Part 11 compliant software package Analysis software	

SPS model	SPS23-100n	SPS11-10 μ
Number of samples & sample size	23 samples \varnothing 18 mm	11 samples \varnothing 50 mm
Temperature range	5 °C to 40 °C (\pm 0.1 K)	5 °C to 60 °C (\pm 0.1 K)
Humidity range*	0 % RH to 98 % RH (\pm 0.5 % RH at 10...30 °C)	
Load range min to max	<3 mg to 2 g	<50 mg to 220 g
Balance resolution	0.1 μ g	10 μ g (up to 111 g) 100 μ g (up to 220 g)
Repeatability**	\pm 1.5 μ g	\pm 20 μ g (up to 111 g) \pm 100 μ g (up to 220 g)
Hardware options		Video camera Raman spectroscopy Water vapor transmission rate Large objects sample tray
Software options	21 CFR Part 11 compliant software package Analysis software	

* Maximum humidity is reached at chamber temperature slightly above room temperature.

** Root Mean Square (RMS). The specified values assume that the system is installed in an environment suitable for the operation of micro- or sub-micro-balances.

Balance stability and repeatability

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Balance calibration

Fully automated calibration procedure with internal reference weights.

Drift compensation

Drift effects from the balance are automatically compensated by an internal procedure using an empty reference pan on the sample carousel.

Balance stability

A drift-free measurement enables excellent repeatability of weighing results independent of the duration of a measurement.

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High sample throughput

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Multi-sample capability

With up to 23 samples measured simultaneously, the ProUmid DVS instruments enable a very high sample throughput compared to a single-sample instrument.

Continuous weighing

By stepwise rotating and lowering of the carousel, the sample pans are placed one after the other on the micro-balance and weighed in predefined time intervals.

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Auswechselbare Probenteller

Exchangeable sample trays

The sample trays for sample pans of different sizes can be changed in just a few seconds.

Variety of samples

All different kinds of samples such as fine powder, granules or tablets, paper and foils, textiles and fibres and also heavy and bulky samples like wood or construction materials can be analyzed with the ProUmid DVS instruments.

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Hardware accessories

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Video camera

Integrated high-resolution CCD camera for detection of moisture induced changes to the samples such as swelling, hygroscopic expansion, liquefaction, as well as changes in color or texture.

Water vapor transmission rate (WVTR)

Permeability Kit for the determination of the water vapor permeability of foils and films according to the EN ISO 7783-1 "dish method for free films".

Raman spectroscopy

Sensor fusion combining water vapor sorption and Raman spectroscopy opens new and exciting perspectives for the analysis of solid materials.

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
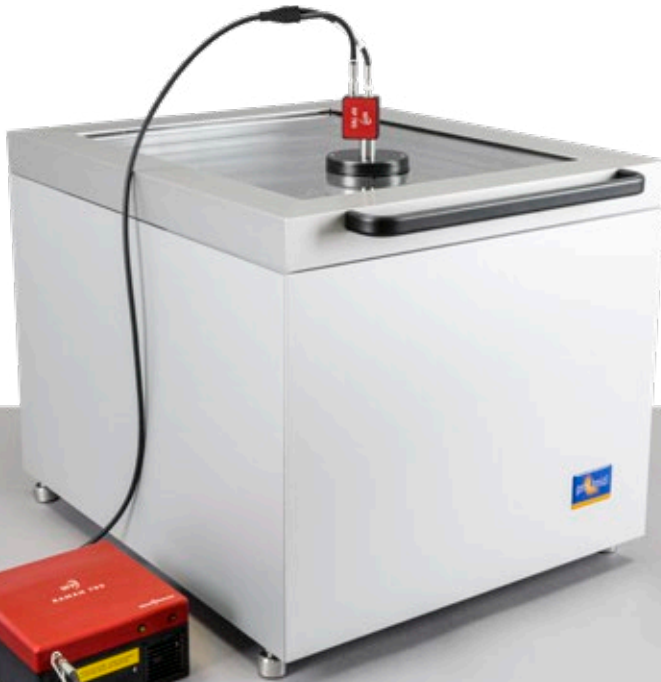
Permeability Kit

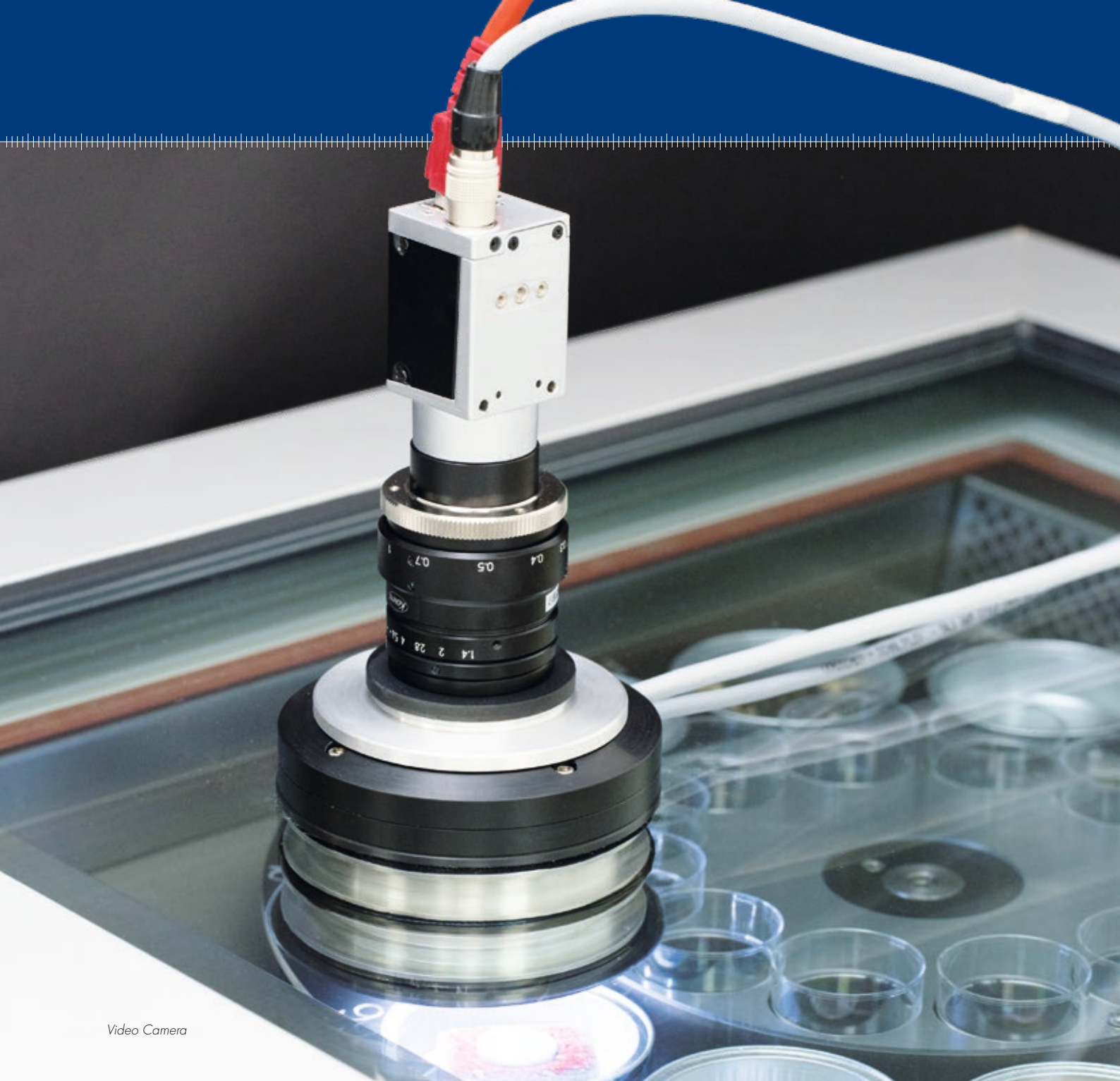Video Camera

Working with the ProUmid DVS instruments

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Sample preparation

Sample material is applied directly into the pans on the sample carousel. The easily accessible measurement chamber enables a convenient sample preparation.

Comparability of results

For perfect comparability of different products in a single measurement, all samples are exposed to the same temperature and relative humidity conditions.

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Fully automated process

Measurement and data acquisition are fully automated. The large glass cover enables the observation of the samples during a measurement. Acquired data can be exported during a running analysis.

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Operating software

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Step-by-step procedure

- User-friendly and intuitive interface
- Programmable sorption/desorption cycles and unlimited number of cycle repetitions
- Humidity cycles as step profiles, ramps or plateaus
- Templates for standard measurements

Data security





Optional software to meet the requirements for data security of electronic records in accordance with 21 CFR part 11.

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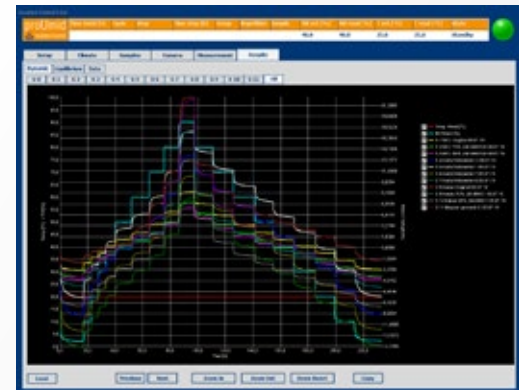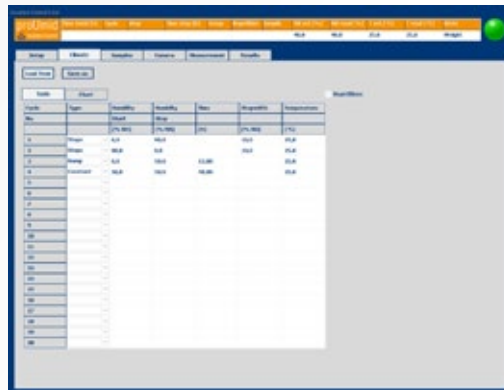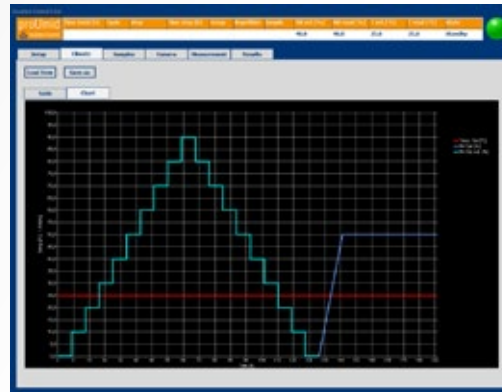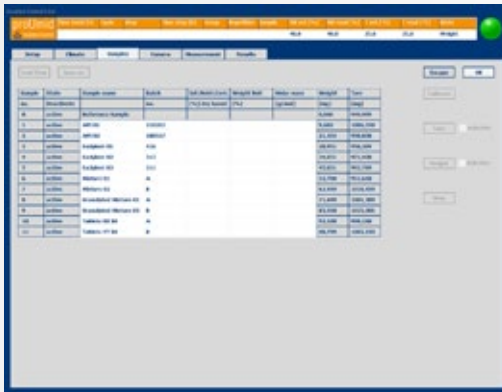Data export and data analysis

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Data export

Acquired data is saved in a binary coded data format. Several options are available for data export into an editable format:

- Import into an Excel® template with fully configurable graphical and tabulated data
- LIMS compatible data format
- Advanced analysis software

Analysis software

Optional software package for data import and advanced data analysis :

- Tabular and graphical display of data
- Configurable and editable diagrams
- Advanced multi-sample data handling
- Curve fitting functionality (BET, GAB, etc.)
- Amorphous content calculation
- Calculation of the water vapor transmission rate
- Calibration tool
- Report generator

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

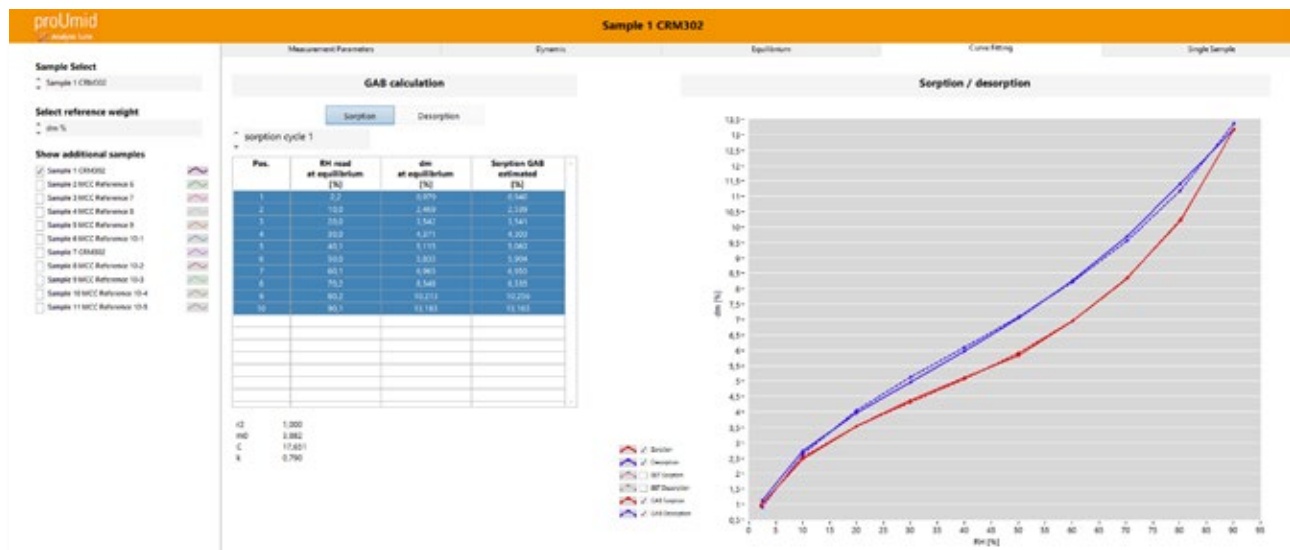
136,617

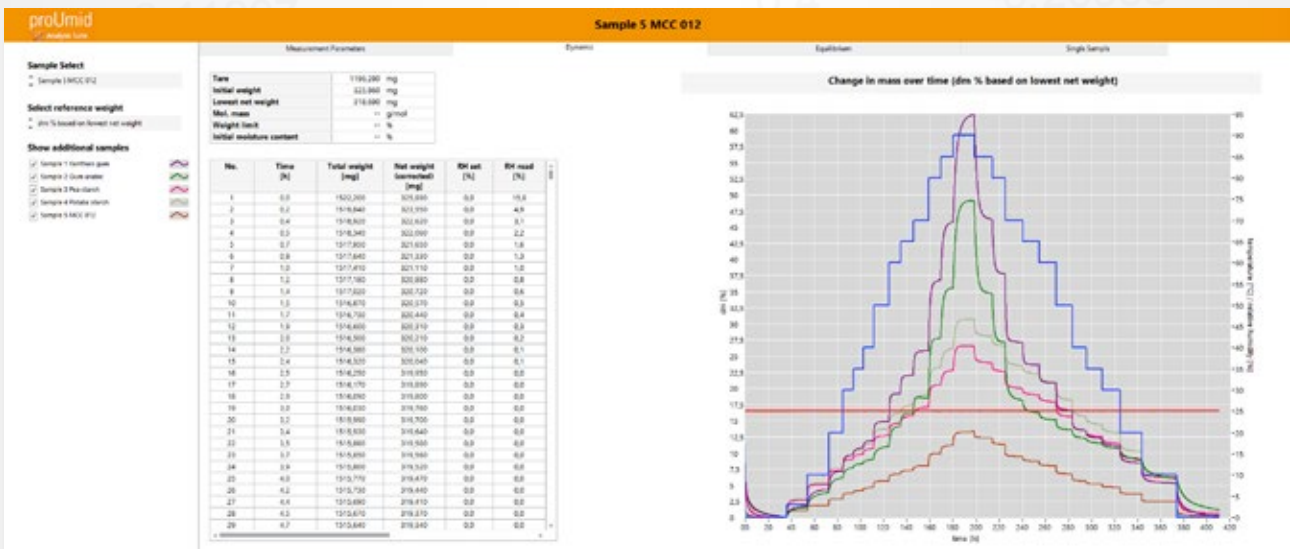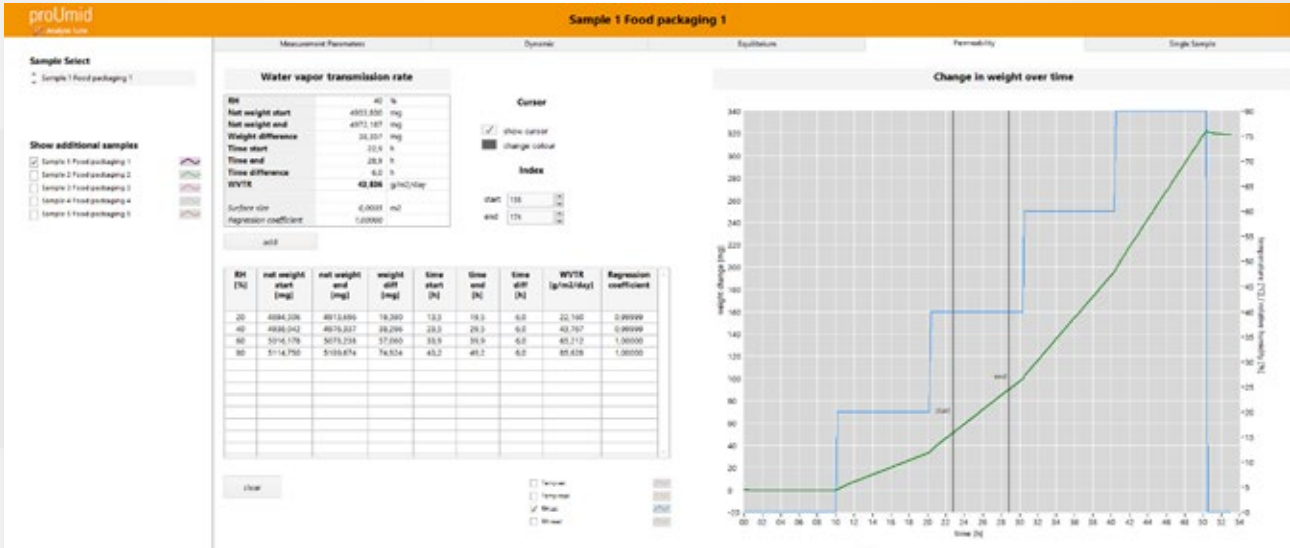Application examples

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DVS is an established standard method in many industries such as pharma, food and feed, fertilizer and agrochemicals, fine chemicals, detergents, cosmetics, building materials, textiles, electronics, archeology.

Dynamic vapor sorption is applied to study:

- Kinetics of moisture sorption/desorption
- Equilibrium water content at different relative humidity levels (sorption isotherm)

Special applications are among others:

- Structural changes (hydrate formation, crystallization)
- Amorphous content determination
- Shelf-life and stability studies
- Deliquescence
- Thermodynamic parameters such as heat of sorption
- Permeability and diffusion studies
- Swelling and hygroscopic expansion
- Quality assurance: determination of minimal differences between batches of a product

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
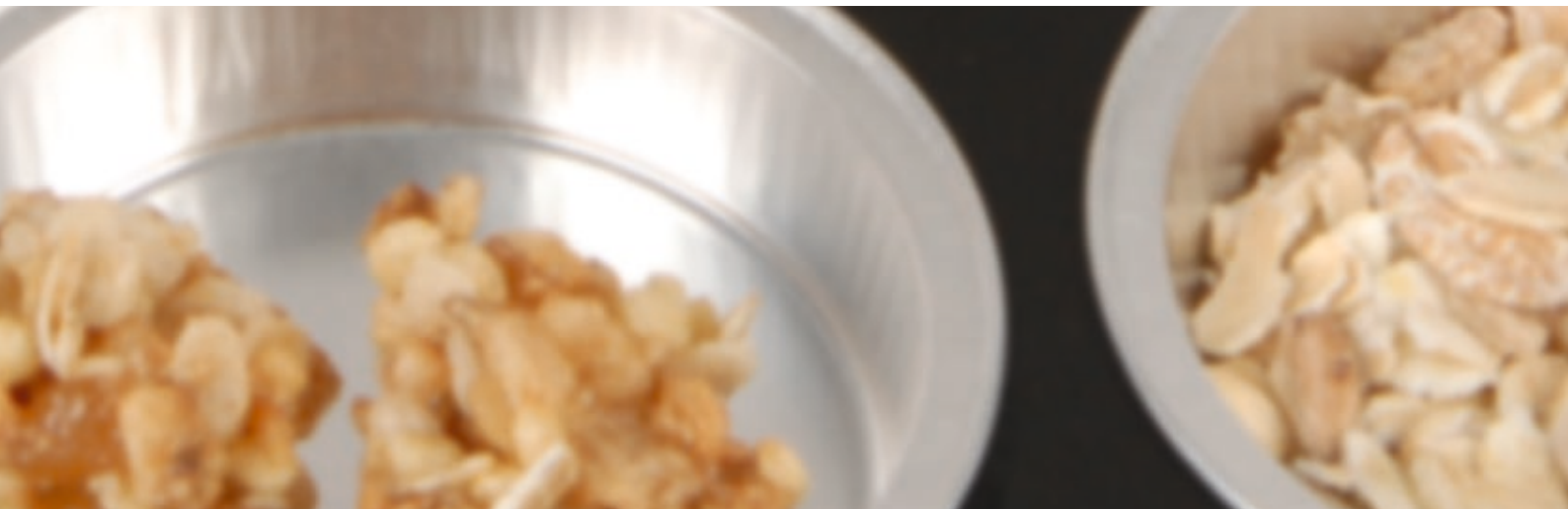
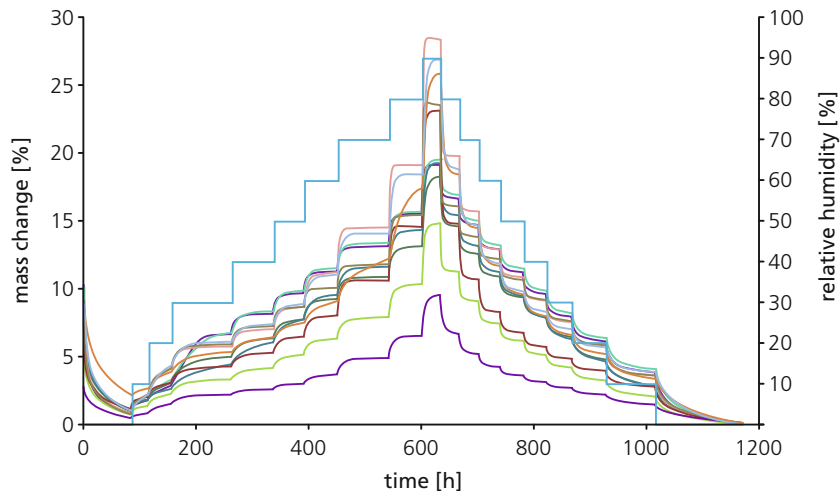
Multi-sample analysis
Longterm study with food materials at 25 °C.

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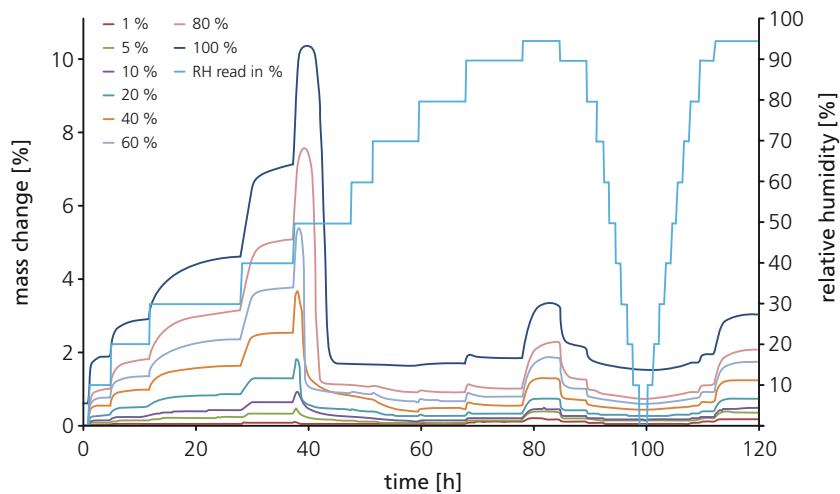Determination of amorphous content

Sorption kinetics of lactose mixtures with different levels of amorphous lactose (crystallization peak at 50 % RH). The data enable to create a calibration curve for amorphous content determination of lactose mixtures.*

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
* Data kindly provided by the University of Innsbruck, Institute for Pharmacy.

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Application examples

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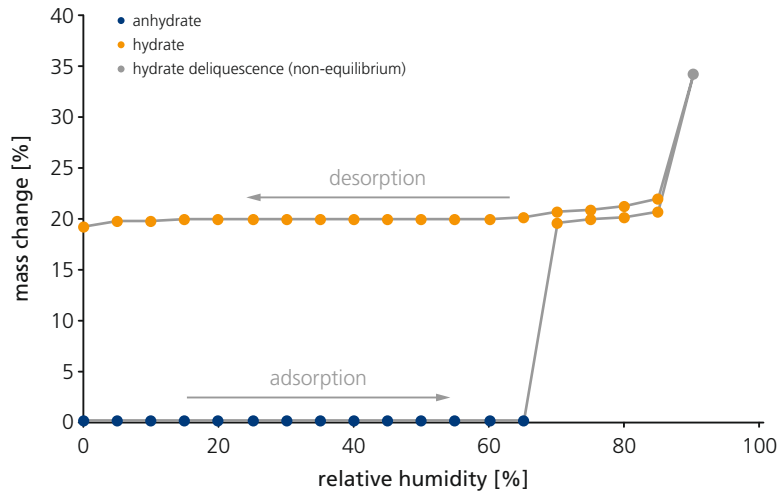
Hydrate formation

Sorption isotherm of L-lysine HCl at 25 °C. Complete adsorption-desorption cycle with transition from anhydrate to hydrate crystals.

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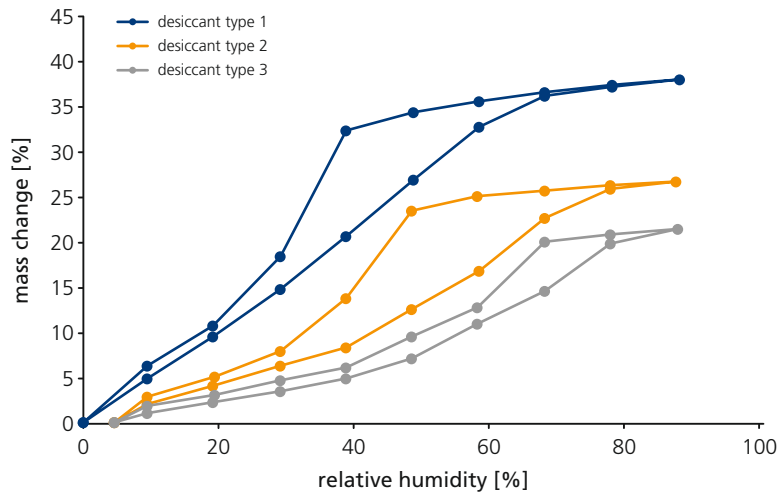
Porosity and Hysteresis

Sorption isotherms of three types of desiccant at 25 °C.

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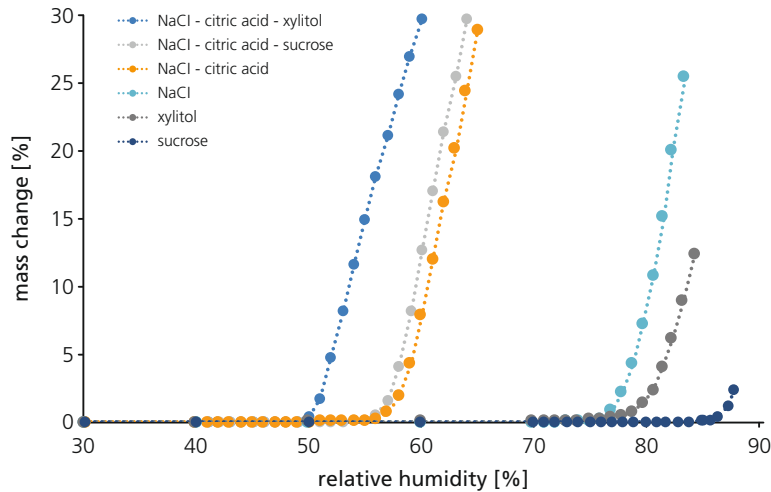
Deliquescence

Deliquescence of water soluble food ingredients and the effect of deliquescence lowering for blends of ingredients.

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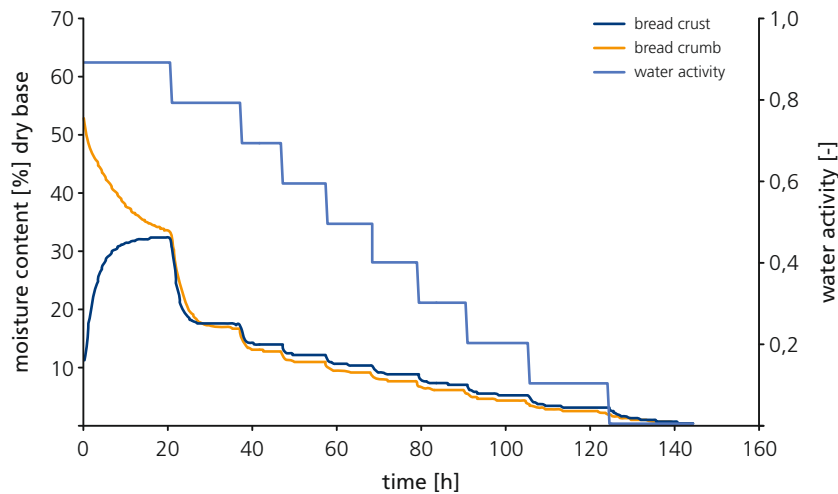Shelf life of food products

The desorption kinetics of fresh bread crust and fresh bread crumb.*

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* Data kindly provided by the Institute of Food, Nutrition and Health, ETH Zürich.

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The SPS in a nutshell

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High sample throughput

Analyzing up to 23 samples simultaneously saves valuable measurement time compared to a single sample instrument and enables the results to be directly compared.

Large dynamic weight range

Different balance resolutions from 0.1 µg for samples of a few mg up to 10 µg for heavy samples enable the analysis of both, tiny samples as well as large samples with lowest water uptake.

Balance stability

The automatic drift compensation enables excellent repeatability of weighing results independent of the duration of a measurement.

Versatility

Optional sensors like Raman spectroscopy or video camera and complementary methods for permeability testing extend the application range of the SPS.

Sustainability

The SPS instruments are designed for a long life cycle. The first device sold in 1999 is still in daily operation and runs with the latest software version.

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ProUmid GmbH & Co. KG